


# Airpulse PK Diagram

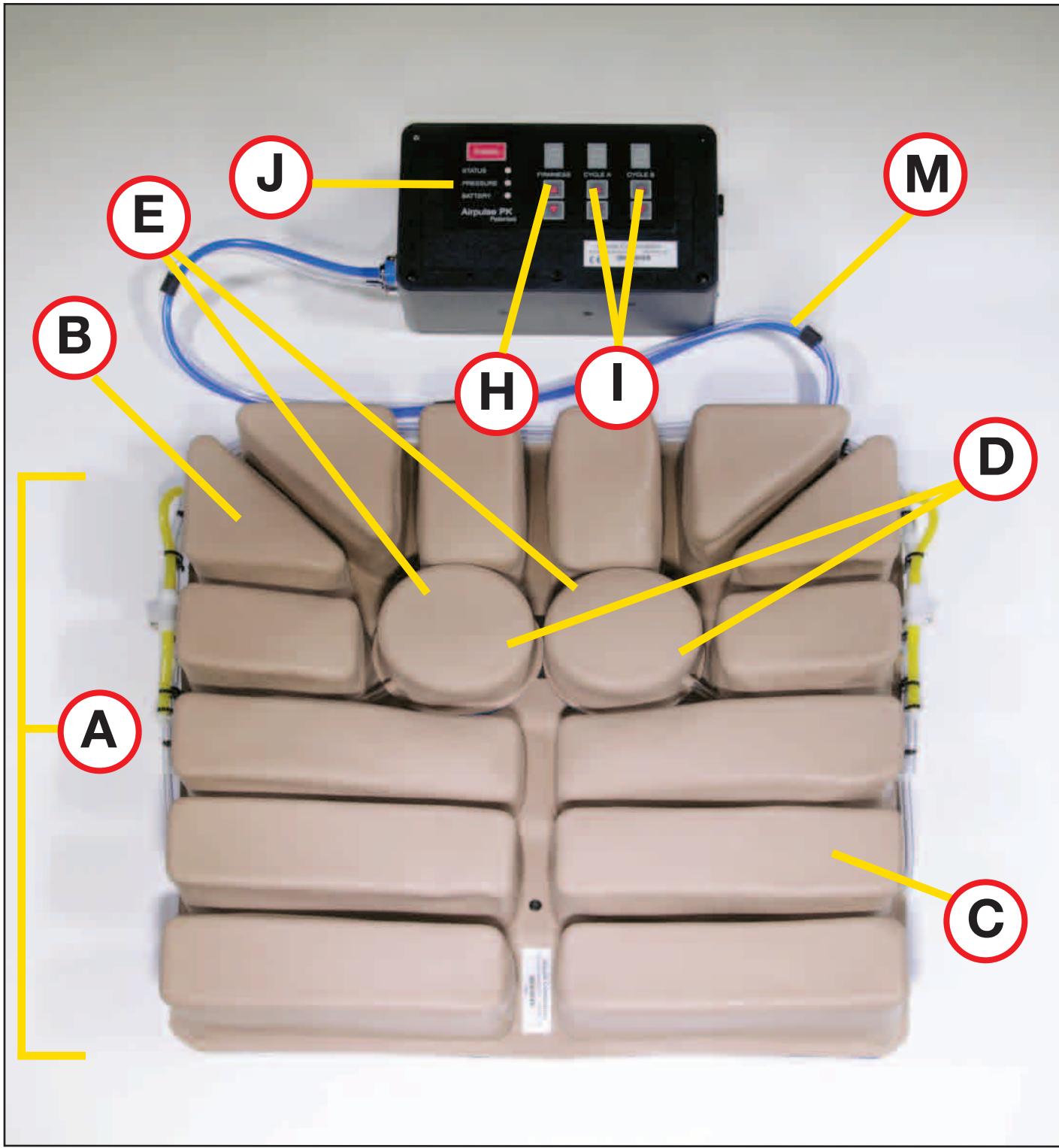

## AQUILA CORPORATION

- A** The cushion is made of soft but durable Neoprene and is custom made to any size including for bariatric situations.
- B** Air cells are sculpted to fit the contour of the body to assist positioning and fit.
- C** Air cells have internal bladders for extra strength and are reinforced with foam for safety to prevent a bottoming out event.
- D** Ischial cells have a Velcro backing and can be repositioned as desired or removed to make room for a wound vac.
- E** The round ischial cells can be easily disconnected so they do not inflate for near zero upward pressure under an existing sore in this critical area.
- F** Covers are included for the cushion, tubes and controller. The cushion cover is made of a stretchable, anti-microbial fabric on the top and non-slip material on the bottom. Incontinent covers are available.
- G** Charger (Double Battery and Hybrid Models) is fully automatic and plugs directly into the controller to charge the battery. No need to remove the battery pack. International models available.
- H** User adjustable firmness setting give you control over cushion inflation pressure.
- I** User adjustable cycle times give you control over frequency of inflation.
- J** Visual and audible indicators for pressure fault and low battery alert.
- K** Power cables (Power chair Model) with a quick disconnect attach directly to the battery terminals and draw power evenly from both chair batteries.
- L** Mini Remote Control (optional) with On/Off switch, system status, low pressure, and battery status indicator lights and beeper.
- M** Air Hose (29 in.) from controller to cushion.
- N** Deluxe Remote Control is an optional accessory that allows you to operate the system without accessing the control box. All system settings can be changed using the Deluxe Remote Control.